

# DID YOU KNOW?

NO: 58

June 7, 2010


## STERILIZER OPTIONS - THE COMPRESSED AIR-OVER COOLING SYSTEM

**Purpose** The Air Over feature provides protection from 'Load Boil Over' when running a cycle with a liquid load

**Parameters** LIQ COOL TIME

**Operation** Air-over cooling incorporates air over-pressure during the cooling phase to ensure that the liquid media does not boil over when steam pressure is removed. Steam is displaced in the chamber with cool air to cool chamber walls, promote cooling of load and shorten exhaust phase of a liquids cycle. Filtered compressed air is supplied by the end user.

**Exceptions** None

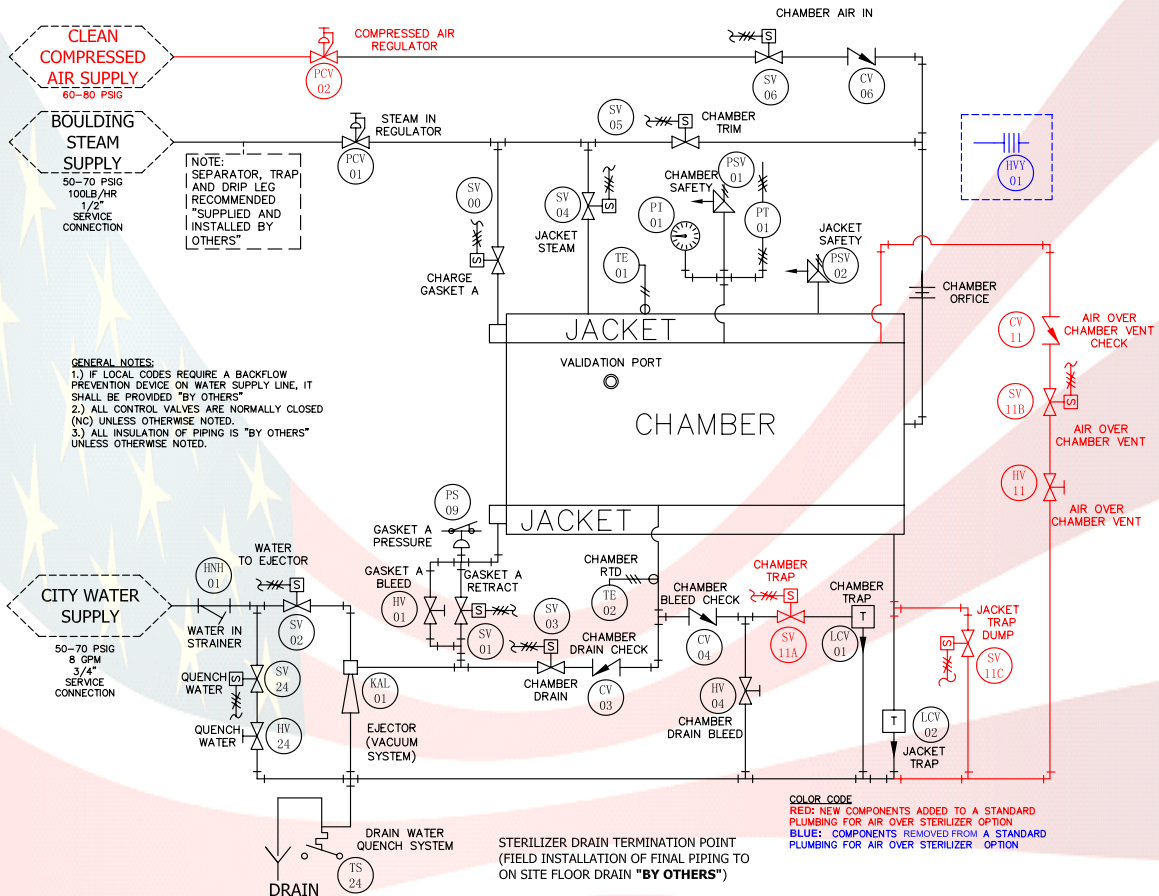

# PRIMUS